


BODYGUARDIAN™ MINI PLUS

Patient Instruction Manual



FULLY WATERPROOF

Featuring a waterproof monitor safe for showering, bathing, and submerging in water as deep as 3 feet.

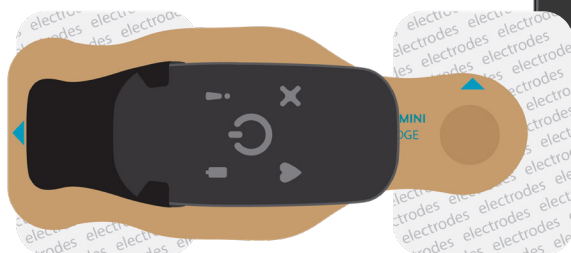

Welcome to **BODYGUARDIAN™ MINI PLUS**

Your physician has prescribed you the **BodyGuardian Mini PLUS** heart monitor and monitoring center service provided by **Boston Scientific Cardiac Diagnostic Services** to detect irregular heart rhythms. You must wear the monitor for the length of time your doctor has prescribed.



Boston Scientific Cardiac Diagnostic Services is not an emergency response service. Contact your doctor immediately if your symptoms worsen. **If you are experiencing a medical emergency, dial 911.**



Did you receive your monitor at home?

If you are receiving the monitor at home and require assistance, call Boston Scientific Cardiac Diagnostics at: **888.500.3522**




When your study concludes

Remove BodyGuardian Mini PLUS as soon as your physician or qualified health care professional ends your service.

Any data or information transmitted or received by Boston Scientific Cardiac Diagnostic Services after your physician or qualified health care professional has ended your service will not be processed, reviewed or analyzed by Boston Scientific Cardiac Diagnostic Services, and Boston Scientific Cardiac Diagnostic Services shall have no liability for any such information or data.

Overview of your monitoring prescription

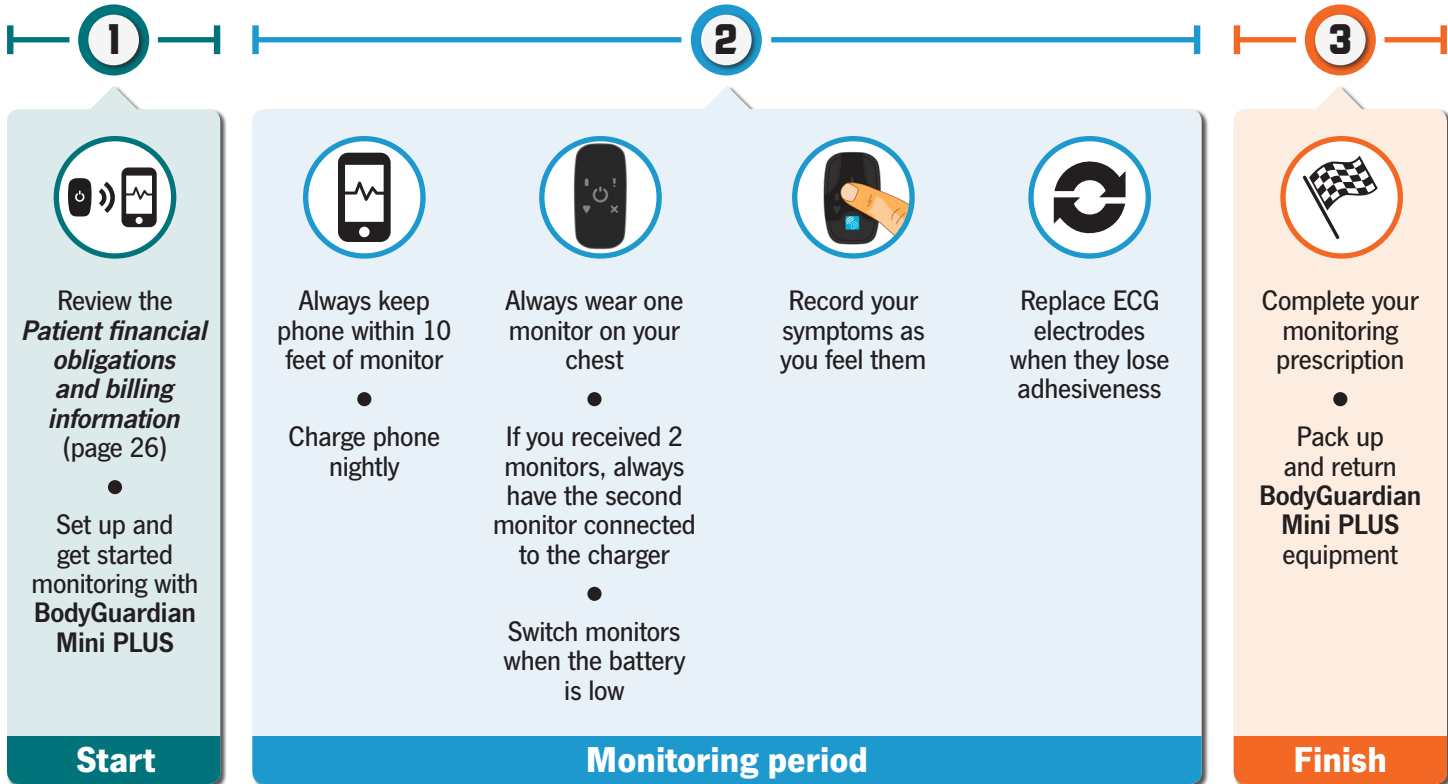

Patient services

- 📞 888.500.3522
- 🕒 24 hours a day, 7 days a week
- ✉️ monitortroubleshooting@CDxBsci.com

Returning your monitor

- 📞 888.500.3522
- 🕒 7:00 a.m. - 7:00 p.m. CST, M-F
- ✉️ monitorrecovery@CDxBsci.com

General information

- Company website: CDx.BostonScientific.com
- Manufacturer's Instructions for Use: manuals.CDx.BostonScientific.com

Billing assistance

- 📞 888.747.4701
- 🕒 8:00 a.m. - 5:00 p.m. CST, M-F
- ✉️ reimbursementinfo@CDxBsci.com

Patient benefit quotes

- 📞 888-747-4760
- 🕒 7:00 a.m. - 5:30 p.m. CST, M-F
- ✉️ patientbenefitquotes@CDxBsci.com



Instructional videos with subtitles to help you use your monitor are available on the BodyGuardian smartphone. On the BodyGuardian phone, tap **Help & Settings** → **Videos**



The BodyGuardian phone cannot be used to make phone calls, including calls to 911, or perform other smartphone functions. Viewing the phone is prohibited while driving, operating heavy machinery or in situations that may jeopardize your safety.


Equipment




BodyGuardian phone



monitors (1 or 2)

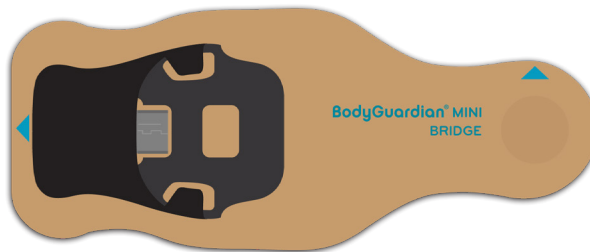

Bridge

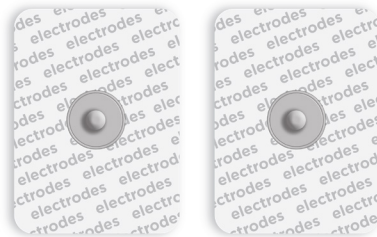

ECG electrodes



Charging cord

(equipment may differ from pictured)



Do not use your own personal device chargers. Only use the charger(s) provided to charge this equipment.

GETTING STARTED

1 Prepare your skin

You can wear BodyGuardian Mini PLUS either vertically or horizontally.



If hair is present, remove the hair where you will place the monitor.




Use soap and water to clean the area.
Dry the area completely.




Do not apply lotions, oils, or perfumes.
Do not use alcohol wipes.
Do not abrade the skin.



Decide which chest placement you will use before you prepare your skin.

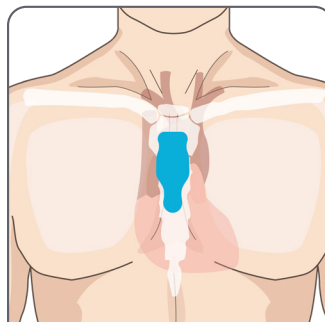


Preferred (vertical)

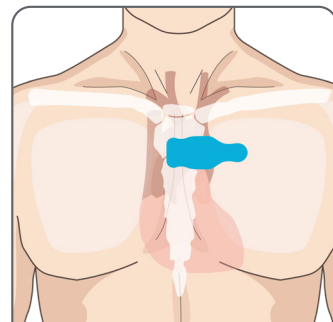

Alternate (horizontal)

2 Snap the ECG electrodes to the Bridge

Press the ECG electrodes firmly onto the Bridge until you hear a snap.




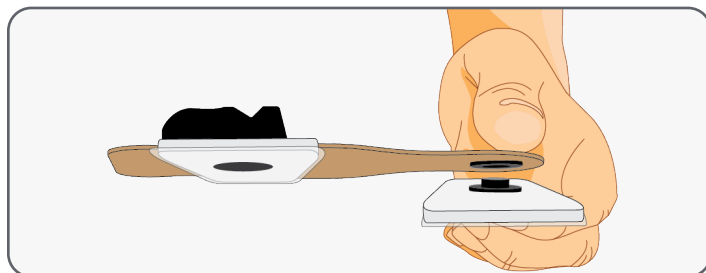


3 Attach the monitor to the Bridge

Slide the monitor into the black molding on the Bridge until the micro-USB connector is inserted completely into the monitor.


Tips

- This is most easily accomplished on a flat surface, like a table.
- To line up the micro-USB connector with the port on the monitor, you may need to wiggle the monitor slightly to align.
- To ensure a good connection, pinch together with thumb and index finger:

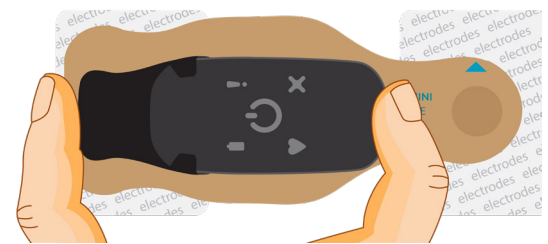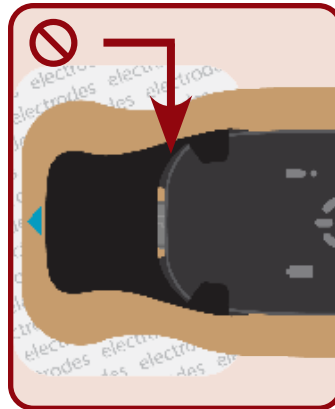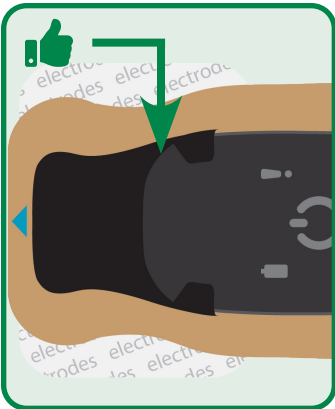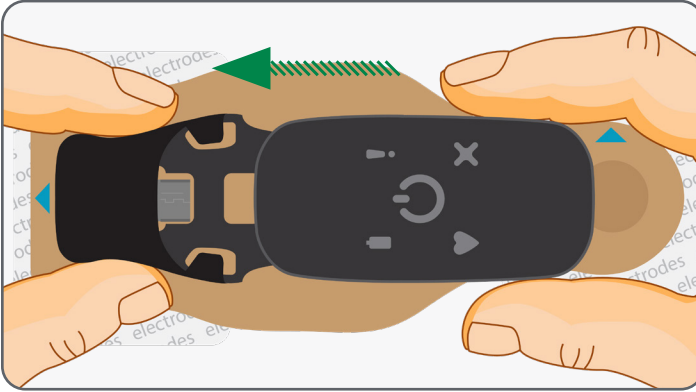


4 Attach to your chest

1. Peel away backing from the electrodes to expose the adhesive.
2. Place on chest in one of the placement locations (**pages 9-10**).
3. Press firmly around the surface of the electrode to ensure adhesion to the skin.



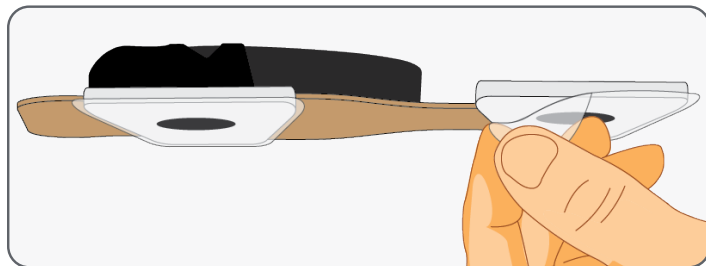


Apply only to intact, clean skin. Do not apply over open wounds, lesions, infected or inflamed areas. If you have skin allergies, hypersensitivities, or experience skin irritation while wearing the monitor, consult your doctor to select the most appropriate option.




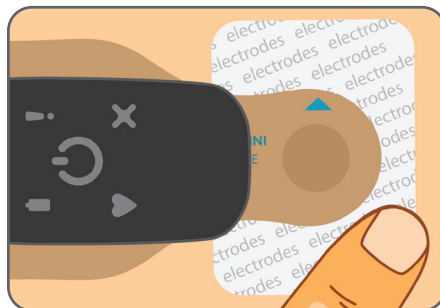


Preferred placement (vertical)

Place the top of the Bridge just below the jugular notch, at the base of the neck at the intersection of the collar bone.

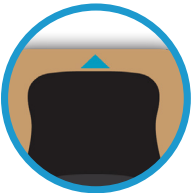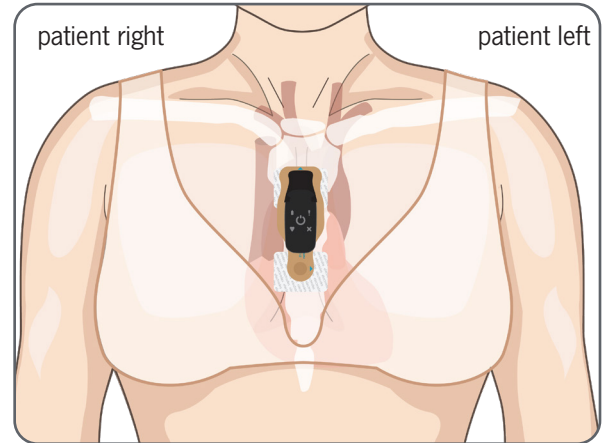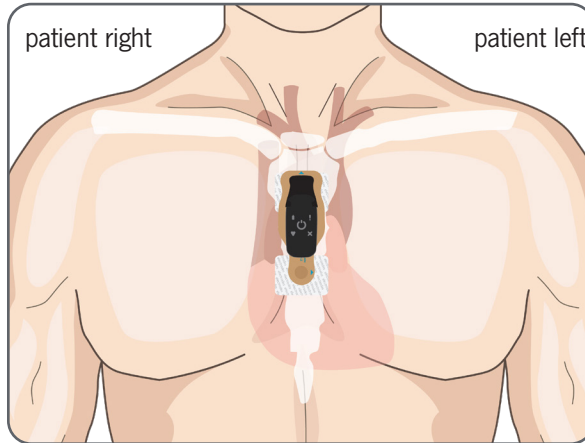



Ensure the arrow at the top of the Bridge on your chest is pointing up.

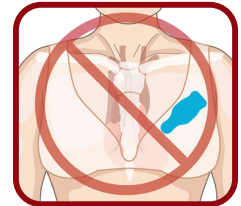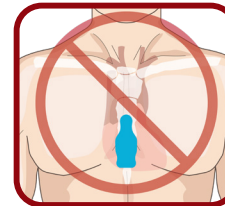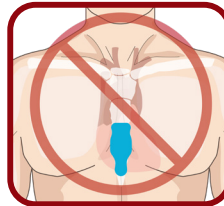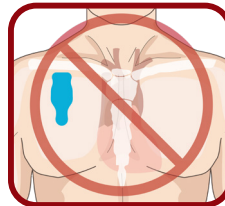







GETTING STARTED

Alternate placement (horizontal)

Place the right side of the Bridge near the center of the chest over the sternum.

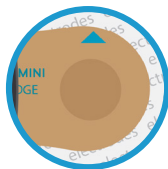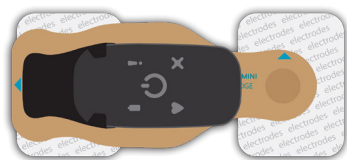

Ensure the arrow on the Bridge over your left chest, near your left arm, is pointing up.

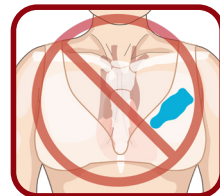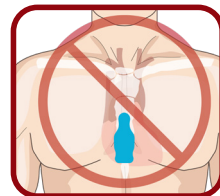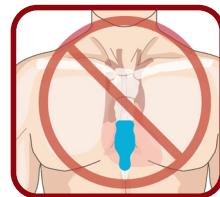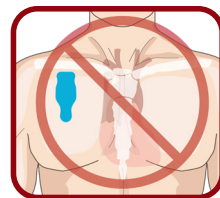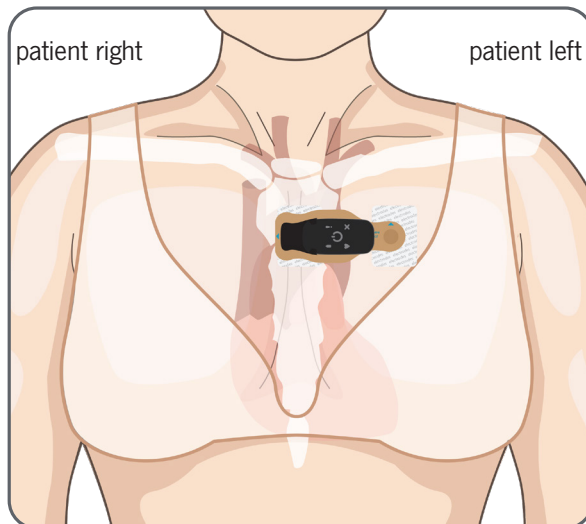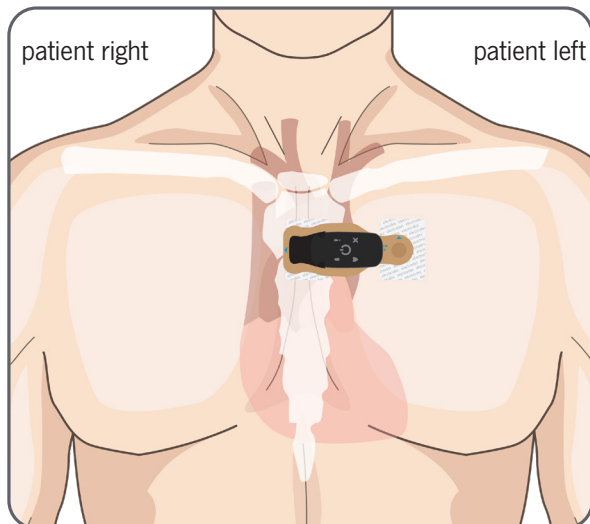


5 Turn on phone

Press and hold the **Power** key located on the right edge or the back of the phone, depending upon the model. Once powered on, the phone displays the

ACTION REQUIRED screen.

6 Turn on monitor


Press and release the center button  on the monitor once. As the monitor successfully powers on, it beeps periodically as it connects to the phone.

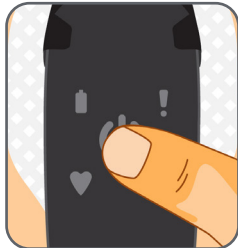

How do I know if the monitor is on?

Visual indicator

An LED light on the monitor blinks every 5 seconds. The monitor is designed to be unobtrusive, so the lights are faint. They can best be detected in a dark or dimly lit room.

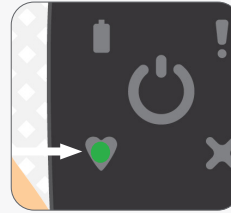


Audible indicator

If you cannot detect the blinking LED light on the monitor, press and release the center button. One beep indicates that the monitor is on. Three beeps indicate that the monitor was off, but is now powering on.


✓ Success

Stay within 3 feet of the phone, remain still, and wait for the phone to display a status of **Monitoring**.

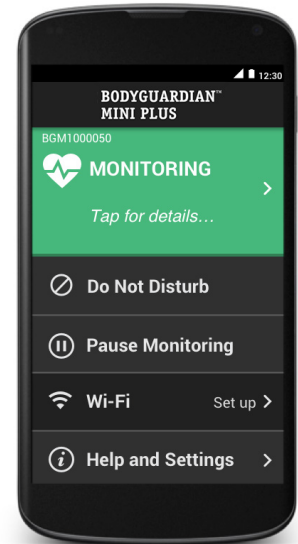

Tip The smartphone supports multiple languages. Tap **Help and Settings** → **Settings** → **Change Language**.

Keep the phone within 10 feet of the monitor on your chest

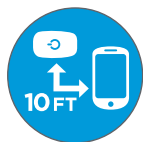



Keep the phone with you at all times. However, if (1) you get separated from the phone or (2) the phone battery depletes, the monitor will continue to collect data and automatically send it to your health care provider when the phone is back in range of the monitor.

If the phone moves out of cellular network range . . .



Your monitoring will continue uninterrupted, but your data won't be sent to your health care provider until you are back in cellular range.

Move to a location with a better cellular signal at your earliest convenience.

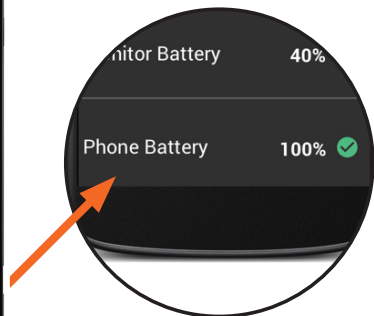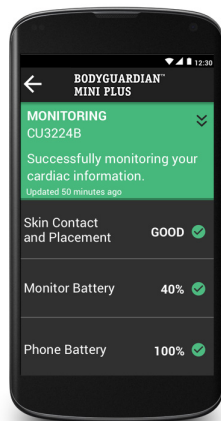

Charge the phone nightly

Ensure the phone is within 10 feet of the location where you sleep. Leave the phone turned on and plug it into the phone charger and into a wall power outlet.

The phone displays its battery status as pictured at left.

Note: If you don't want to be disturbed by alerts while sleeping or while at an event, tap **Do Not Disturb** on the phone.

Wi-Fi

You can connect the phone to an available Wi-Fi network. This can be particularly useful if you find yourself in an area of poor cellular coverage.



1. On the home screen, tap **Wi-Fi**.
2. Ensure the Wi-Fi button is toggled to **ON**.
3. From the list displayed, tap a network that you trust, then enter the password (if required).

Record your symptoms


If you begin to feel symptoms related to the reason your monitor was prescribed such as **dizziness**, **chest pain** or **shortness of breath**, you can manually record these events.

1. Press and release the center button on the monitor.



After pressing the button, it can take up to 60 seconds for the symptoms prompt to appear on the phone.


2. Hold as still as possible for 60 seconds.
3. On the phone, select the symptoms you are feeling. You may select up to 3 symptoms.
4. Tap **Save**.

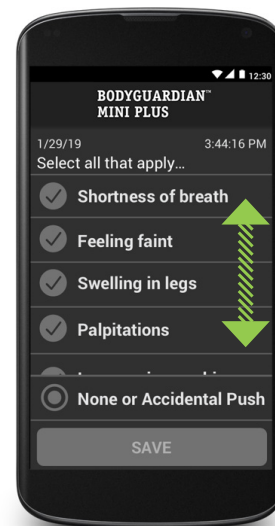

You must scroll to see the full list of symptoms.

When your monitor battery is low

When the monitor battery is low, the phone will alert you.

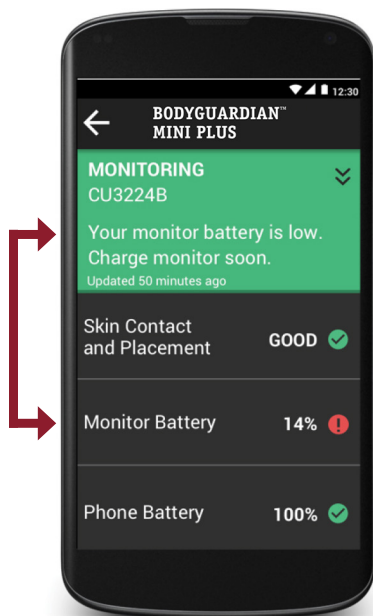


If you have 2 monitors . . .

If you have 2 monitors, switch monitors. Always wear one of the monitors while the other monitor remains connected to the charger.

1. Detach the monitor from the Bridge on your chest.
You do not need to remove the electrodes from your chest to switch monitors.
Do not pause or power off the monitor. The monitor will automatically power off when it is done transferring data and charging.

2. Firmly plug the charger into the monitor.

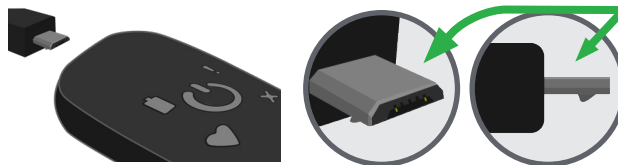

Important

Ensure the micro-USB connector is oriented this way when inserting into the monitor.

! You must plug the charger into the monitor you just removed from your chest before you attach the other, fully-charged monitor onto your chest.


⊘ *Do not use your own personal device chargers. Only use the charger(s) provided.*

3. Attach the charged monitor to the Bridge, then re-attach the Bridge to the electrodes on your chest.
4. Press and release the center button on the monitor to power it on.
 - **Do not hold the button down.**
 - Remain near the phone when attaching the fully-charged monitor.
 - The phone will automatically detect the new monitor.

If you have 1 monitor . . .

If you have 1 monitor, remove and charge the monitor.

1. Remove the Bridge and electrodes from your chest, then remove the monitor from the Bridge.
2. Firmly plug the charger into the monitor.

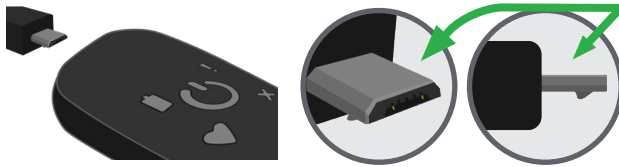

Important

Ensure the micro-USB connector is oriented this way when inserting into the monitor.

 Do not use your own personal device chargers.
Only use the charger(s) provided.

-  **AS SOON AS THE MONITOR IS FULLY CHARGED (1.5-2.5 HOURS),** use *new* electrodes to re-apply BodyGuardian Mini PLUS to your chest.
- See **pages 6-11** for instructions.

Charging the monitor



While charging, the blue light on the monitor blinks.



Once fully charged, the blue light remains solidly lit.



Charging time: 1.5-2.5 hours

Replacing the ECG electrodes



Change the ECG electrodes every 3 days OR when they no longer adhere to the skin. Electrodes are disposable.

1. Remove the monitor and Bridge from your chest.
2. Slowly peel back and roll the adhesive on the back of the electrode away from your chest.
Tip: Dampen with water to ease removal from your skin or remove while in the shower.
3. Re-attach BodyGuardian Mini PLUS to your chest (see **pages 8-10**).

! Do not rip or tear off the electrodes quickly. Improper removal may cause skin irritation.

Water resistance



The accessories (Bridge and electrodes) are waterproof and can safely be submerged in water as deep as 3 feet. The monitor is waterproof when plugged into a body-worn accessory that does not leave the micro-USB port exposed. The micro-USB port should **not** be exposed to water.



The phone is NOT waterproof. Keep away from water.

- You do not need to remove BodyGuardian Mini PLUS when showering or bathing.
- If you will be swimming in water deeper than 3 feet, remove the monitor.
- If micro-USB connector on an unattached monitor is exposed to water, open air dry for at least 1 hour before reattaching to the Bridge.



After showering, bathing, swimming while wearing BodyGuardian Mini PLUS:

1. Pat the Bridge with a towel to dry it off.
2. Press and hold the electrodes against your chest to ensure good contact with your skin.

MRI



BodyGuardian Mini PLUS cannot be worn while undergoing an MRI.

1. On the phone, tap **Pause Monitoring**.
2. Remove **Bridge, monitor** and **electrodes** from chest.
3. After the MRI, snap new electrodes to the Bridge and re-attach to your chest (see **pages 8-10**).
4. Tap **Resume** on the phone.

Air travel



BodyGuardian Mini PLUS is safe to wear while traveling on aircraft. You do not need to remove the monitor. If you will be traveling on aircraft:

1. Contact Boston Scientific Cardiac Diagnostic Services prior to departure to document dates of travel.
2. Power off phone after boarding.
3. Power on phone after landing.

Monitor indicators

This section explains the sounds and LED lights on the monitor. The phone communicates this same information and is the best method for understanding the status of your monitoring session.

If you are curious about what the monitors sounds and lights mean, you can use this information as a reference.

Note The monitor lights are faint. They can best be detected in a dark or dimly lit room.




Errors				
Sounds	Lights*		What it means	Caused by
♪♪	and	✖ every 1 sec	Poor skin contact.	Monitor cannot detect good contact with skin to begin monitoring.
none	and	✖ every 5 sec	Poor skin contact.	Monitor cannot detect good contact with skin during already established monitoring session. If this persists for 30 minutes or more, replace your electrodes.
♪♪	and	✖ every 5 sec	Poor skin contact.	Monitor has not detected good contact with skin for over 8 hours. Replace your electrodes.
none	🔋!		Charging issue.	Battery charge not reaching target level in pre-defined time.

Warnings

Sounds	Lights*	What it means	Caused by
 and 	every 5 sec	Monitor battery is low.	Low monitor battery. Charge the monitor (page 15).

Informational

Sounds	Lights*	What it means	Caused by
none	 every 5 sec	Monitoring successfully.	
none	 every 5 sec	Monitoring successfully.	Monitor is communicating with phone over Bluetooth.
 and none		Symptom marked.	Press-and-release button push.
 and 		Monitor is powering on.	Press-and-release button push.
 and    		Monitor is powering off.	Press-and-hold (10-12 sec) button push.
After 10-12 seconds, all four monitor lights flash once, followed by a long beep.			

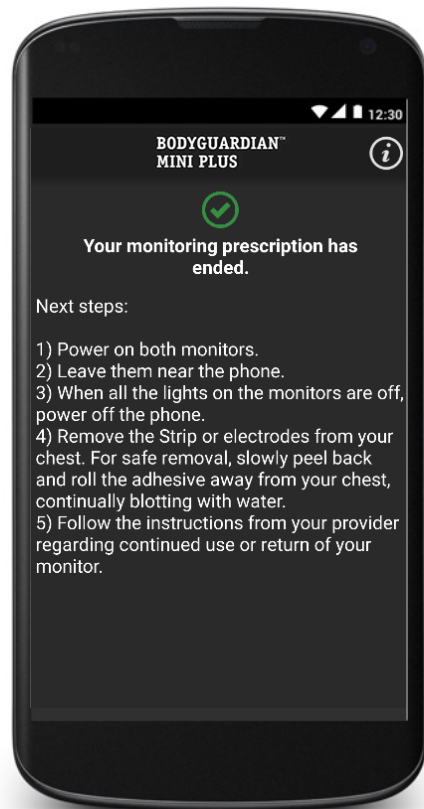

When your prescription ends

The phone will alert you when your prescription has ended.

When you see this message on the phone:

1. Remove the monitor on your chest, then place it near the phone.



Do not power off the monitor.

2. Power on the other monitor and leave it near the phone.
3. Wait for the monitor(s) to power off. The monitors will power off automatically once they are done sending data. You will know the monitor is off when none of the four lights are lit.



4. Power off the phone.
5. Pack up and return your equipment (see **page 21**).

Return your equipment immediately after your service ends

If you are returning your equipment to your health care provider:


1. Follow the instructions provided by your health care provider.
-


If you are returning your equipment to Boston Scientific Cardiac Diagnostic Services:

1. Place the powered-off monitor(s) into the foam pouch provided, then into the original postage-paid box. Then place the phone, charging cord, any unopened Bridges and ECG electrodes, and patient survey into the box.
2. Remove the adhesive strip from flap on the shipping box and seal securely. **The prepaid return label is already affixed to the bottom of the box.**
3. Drop package off at the shipping carrier listed on the prepaid return label. For assistance, call Boston Scientific Cardiac Diagnostics at **888.500.3522**.



If you have lost or damaged the equipment, Boston Scientific Cardiac Diagnostics will mail you an invoice for the cost.

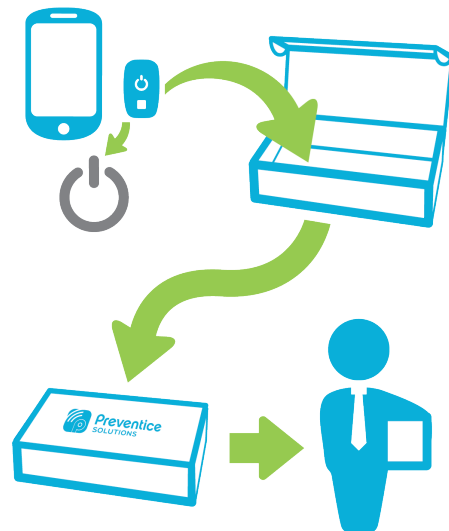

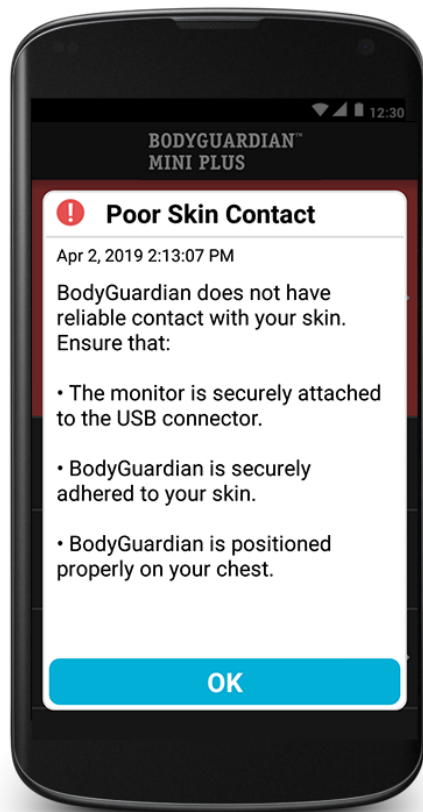

***Poor Skin Contact* message appears on the phone**

CAUSE

This message indicates poor contact between your skin, the monitor and Bridge (or electrodes). It is NOT related to the Bluetooth or cellular connection of the phone.

RESOLUTIONS

1 Tap OK

Tap **OK** on the phone to dismiss the message.

2 Check the monitor attachment

Ensure the monitor is securely attached to the micro-USB connector.

3 Check the chest application

a. Adhesion

Ensure that the electrodes are securely adhered to the skin without wrinkles or bubbles. When applying the electrodes, apply pressure to ensure good contact with the skin.

b. Chest placement location

Ensure that BodyGuardian Mini PLUS is positioned properly on the chest. See **pages 9-10** for chest placement locations.

Tip: For the **horizontal Bridge position**, some patients acquire a better signal when the monitor is placed higher up on the chest, closer to the collarbone, to avoid fleshy tissue areas. BodyGuardian receives a higher quality signal when placed over firmer areas.

c. Skin preparation

Ensure that the skin on the placement location on your chest has been prepared appropriately:


- Remove hair, excess oils, lotions, perfumes, etc.
- Use soap and water to clean the area. Do not use alcohol wipes. Dry the area completely.
- Do not abrade the skin.

4 Replace the ECG electrodes

1. Ensure that the chest placement location is properly prepared. See **3c. Skin preparation**.
2. See **page 16** for instructions on replacing the ECG electrodes.



If the message persists, contact Boston Scientific Cardiac Diagnostics: **888.500.3522**

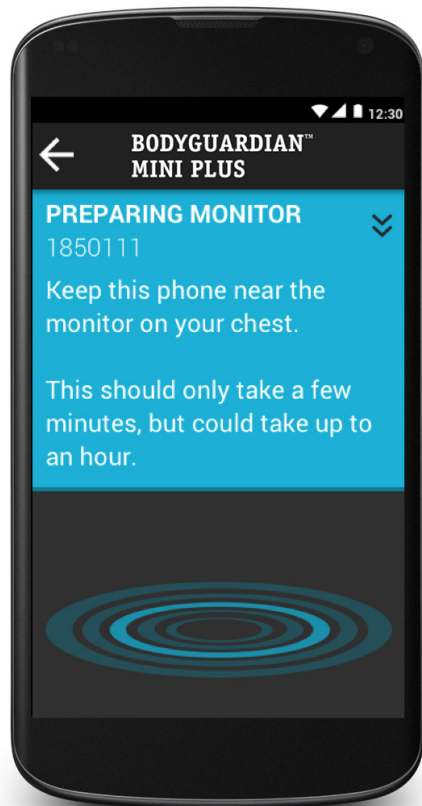

When switching monitors, the phone displays ***PREPARING MONITOR*** for an extended time

CAUSE

This occurs when there is a large amount of collected data on the monitor that is currently being downloaded to the phone.

RESOLUTION

Wait for the phone to finish downloading.

1. Keep the phone within 3 feet of both monitors to allow the downloading process to complete.
2. Once the phone has finished downloading, the status on the phone will change to ***Monitoring***. This may only take a few minutes, but could take up to an hour, depending on the amount of data to be downloaded from the monitor to the phone.



- or -



The phone touch screen dims or shuts off automatically

CAUSE

Even when the phone is collecting data, the phone turns off the touch screen to conserve battery when you do not use the phone for a specified period. This is normal. The phone is still performing its monitoring functions. It has simply turned off the touch screen display.

RESOLUTION

No action is necessary. If you would like to wake the phone screen:

1. Press the **Power** key located on the right edge or the back of the phone, depending upon the model.



A number pad prompting for a passcode appears on the phone

CAUSE

The phone may display a security passcode screen.

RESOLUTION

This screen can be dismissed.

Tap the white arrow in the upper-left corner of the screen.

Patient financial obligations and billing information



On behalf of your physician, you are being provided with a heart monitor from Boston Scientific Cardiac Diagnostic Services. This test will be billed in two parts:

1. Your medical professional will bill your insurance for the in-office hook-up and the final reading of the test results.
2. Boston Scientific Cardiac Diagnostic Services will bill your insurance for the use of the monitor along with providing the physician with 24/7 monitoring center service and all requested data and reports.



Your insurance company will usually send you an Explanation Of Benefits (EOB) describing the amount paid and the amount you owe Boston Scientific Cardiac Diagnostic Services. An EOB is not a bill.

Our goal is to help your physician diagnose your heart condition without delay or concerns about billing. We offer patient-friendly payment options, including a discounted rate at the time of service or a reduced rate when your insurance indicates a large fee payable by you for our services or the physician's services.



All equipment is the sole property of Boston Scientific Cardiac Diagnostic Services. To avoid financial liability and to ensure that you do not receive a bill for the value of the equipment in your possession (up to \$2,400), please return the equipment immediately after your service ends. Any request for an extension of your study beyond your prescribed study needs to come from your physician through a new order sent to Boston Scientific Cardiac Diagnostic Services. Be aware that the second study may not be covered by all insurance providers and you may be responsible for those charges.

Our Notice of Privacy Practices is available for your review at: **CDx.BostonScientific.com**

PATIENT SURVEY

We appreciate your feedback. Please return this survey inside the box once your study is completed. This survey can also be completed online at: **CDx.BostonScientific.com/us/en/patients/patient-support/patient-survey.html**

Patient name _____ City _____ State _____

Doctor's office / hospital _____ Contact me via: ☐ Phone _____

My hook-up took place at: ☐ Home ☐ Physician's office ☐ Email _____

BODYGUARDIAN™ MINI PLUS

	Excellent	Average	Poor	Not applicable
Rate our service level	1	2	3	4
General satisfaction with Boston Scientific Cardiac Diagnostics	☐	☐	☐	☐
Ability to contact Boston Scientific Cardiac Diagnostics at 888.500.3522	☐	☐	☐	☐
Attention given to your comments	☐	☐	☐	☐
Knowledge, professionalism and courtesy of representatives	☐	☐	☐	☐
Explanations given to hook up your monitor	☐	☐	☐	☐
Amount of time the representative spent with you	☐	☐	☐	☐
Rate the contents of the materials	1	2	3	4
Printed instruction manual included in your box	☐	☐	☐	☐

continue ➔



PATIENT SURVEY

BODYGUARDIAN™ MINI PLUS

	Excellent	Average	Poor	Not applicable
Rate your experience using the monitor and its accessories	1	2	3	4
Setting up the monitor	☐	☐	☐	☐
Wearing and changing the monitor	☐	☐	☐	☐
Charging the battery	☐	☐	☐	☐
Adherence and comfort of the Bridge and electrodes	☐	☐	☐	☐
Shipment of supplies for the length of your study	☐	☐	☐	☐

Comments

☐ I am interested in potentially sharing my patient experience and consent to Boston Scientific Cardiac Diagnostics, Inc. contacting me for such purposes at the contact information I provided above, which may include telephone or email.

