

BODYGUARDIAN™ MINI

Patient Instruction Manual



FULLY WATERPROOF

Featuring a waterproof monitor safe for showering, bathing, and submerging in water as deep as 3 feet.



Welcome to **BODYGUARDIAN™ MINI**



Your physician has prescribed you the **BodyGuardian Mini** heart monitor and monitoring center service provided by **Boston Scientific Cardiac Diagnostic Services** to detect irregular heart rhythms. You must wear the monitor for the length of time your doctor has prescribed.

When your study concludes



Remove BodyGuardian Mini as soon as your physician or qualified health care professional ends your service.




Important

Once your prescription has ended, turn off the monitor by pressing and holding the center button  for 9-12 seconds. After 9-12 seconds, the monitor lights will flash, followed by a long beep. This indicates that the monitor has powered off.

Any data or information transmitted or received by Boston Scientific Cardiac Diagnostic Services after your physician or qualified health care professional has ended your service will not be processed, reviewed or analyzed by Boston Scientific Cardiac Diagnostic Services, and Boston Scientific Cardiac Diagnostic Services shall have no liability for any such information or data.

Rx ONLY Caution: Federal law (USA) restricts this device to sale by or on the order of a physician.

Overview of your monitoring prescription

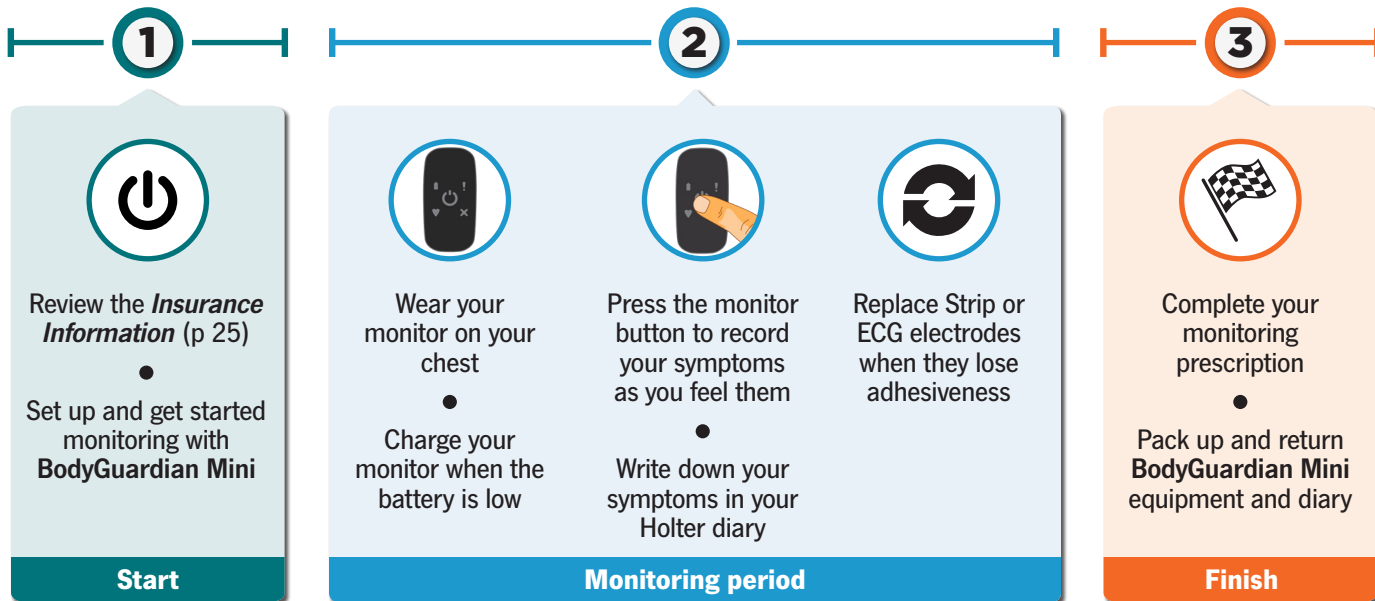

CONTACT INFORMATION

Patient services

- 📞 888.500.3522
- 🕒 24 hours a day, 7 days a week
- ✉️ monitortroubleshooting@CDxBsci.com

Returning your monitor

- 📞 888.500.3522
- 🕒 7:00 a.m. - 7:00 p.m. CST, M-F
- ✉️ monitorrecovery@CDxBsci.com

General information

Company website: CDx.BostonScientific.com
Manufacturer's Instructions for Use: manuals.CDx.BostonScientific.com

Billing assistance

- 📞 888.747.4701
- 🕒 8:00 a.m. - 5:00 p.m. CST, M-F
- ✉️ reimbursementinfo@CDxBsci.com

Patient benefit quotes

- 📞 888-747-4760
- 🕒 7:00 a.m. - 5:30 p.m. CST, M-F
- ✉️ patientbenefitquotes@CDxBsci.com



Boston Scientific Cardiac Diagnostic Services is not an emergency response service. Contact your doctor immediately if your symptoms worsen. **If you are experiencing a medical emergency, dial 911.**



monitor

monitor charging
cord

Do not use your own personal device chargers. Only use the charger(s) provided to charge this equipment.

Your equipment will include either  **Strips** or an  **ECG lead set (with electrodes)**. You will not receive both. Follow the instructions for the equipment you receive.

Strips**ECG lead set****ECG lead wires**

2-wire

or



4-wire

ECG electrodes

Saline wipes
(included with
Strips only)

When you feel a symptom:

1. On the monitor, press and release the **OK** button **O** once. Do not hold the button **O**.
2. Record the **DATE** and **TIME**.
3. Select at least one **SYMPTOM**.
4. Select your **ACTIVITY** at the time you felt the symptom.

	DATE	TIME	SYMPTOM
1	__/__/__	__:__:__ AM/PM	☐
2	__/__/__	__:__:__ AM/PM	☐
3	__/__/__	__:__:__ AM/PM	☐
4	__/__/__	__:__:__ AM/PM	☐
5	__/__/__	__:__:__ AM/PM	☐
6	__/__/__	__:__:__ AM/PM	☐
7	__/__/__	__:__:__ AM/PM	☐
8	__/__/__	__:__:__ AM/PM	☐

Holter diary

Did you receive Strips or an ECG lead set?

This section (pages 6-11) includes instructions for using  **Strips** to wear your BodyGuardian Mini monitor. If you received an  **ECG lead set**, see page 12.




Strip




ECG lead set

 Decide which chest placement you will use before you prepare your skin.

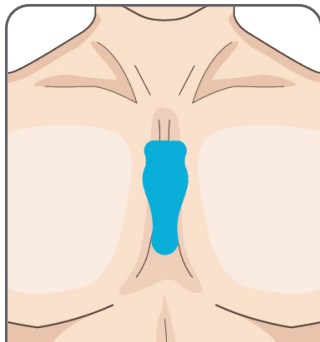

Preferred (vertical)

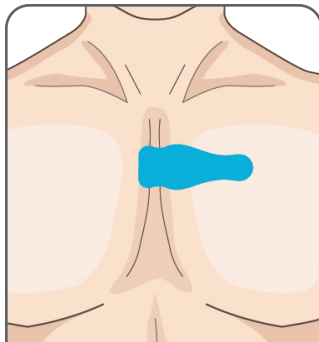

Alternate (horizontal)

1 Prepare your skin

You can wear BodyGuardian Mini either vertically or horizontally.




If hair is present, remove the hair where you will place the monitor.



Use an enclosed saline wipe or soap and water to clean the area. Dry the area completely.



Do not apply lotions, oils, or perfumes.
Do not use alcohol wipes.
Do not abrade the skin.

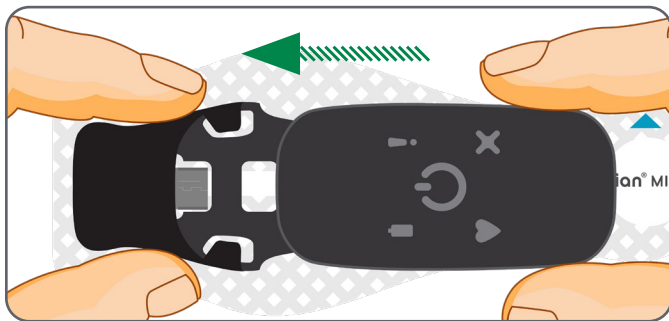




2




Attach the monitor to the Strip

Slide the monitor into the black molding on the Strip until the micro-USB connector on the Strip is inserted completely into the monitor.

Tips

- This is most easily accomplished on a flat surface, like a table.
- To line up the micro-USB connector with the port on the monitor, you may need to wiggle the monitor slightly to align.
- To ensure a good connection, pinch together with thumb and index finger:

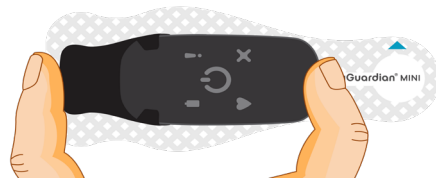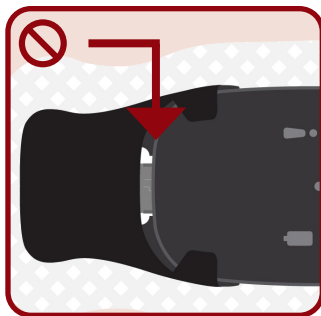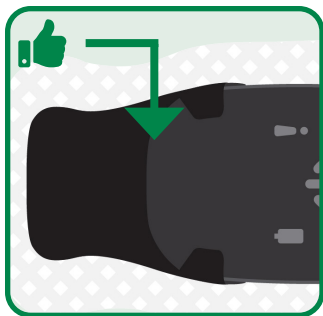

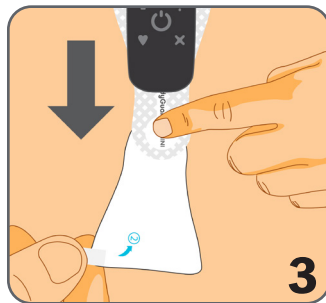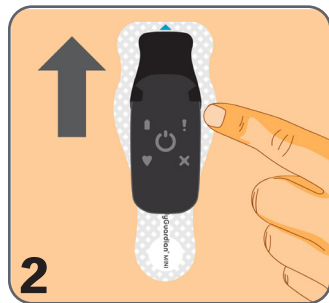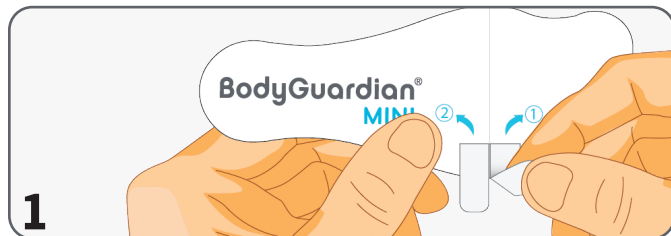

3

Attach to your chest

1. Peel away backing ① from the Strip.
2. Place the portion of the Strip with the adhesive exposed flat on your chest in one of the placement locations (p. 9-10).
3. Slowly peel away backing ② from the Strip as you press the Strip flat on your chest. Press the entire surface of the Strip against the skin to ensure it adheres.



Apply only to intact, clean skin. Do not apply over open wounds, lesions, infected, or inflamed areas. If you have skin allergies, hypersensitivities, or experience skin irritation while wearing the monitor, consult your doctor to select the most appropriate option.




Preferred placement (vertical)

Place the top of the Strip just below the jugular notch, which is the bone at the base of the neck at the intersection of the collar bone.



Ensure the uppermost arrow on the Strip on your chest is pointing up.

❗ Once the Strip makes contact with your skin, you cannot adjust its position.

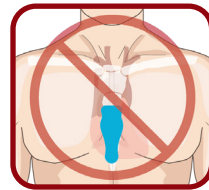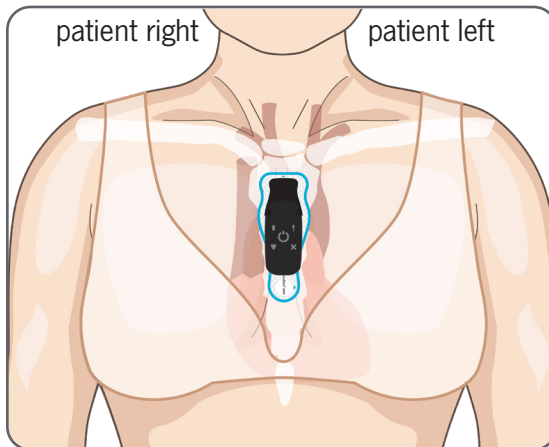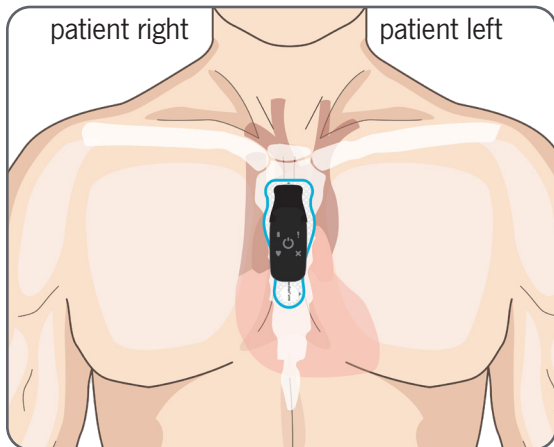

continue →





Alternate placement (horizontal)

Place the right side of the Strip near the center of the chest over the sternum.



Ensure the arrow on the Strip over your left chest, near your left arm, is pointing up.

❗ Once the Strip makes contact with your skin, you cannot adjust its position.

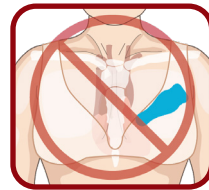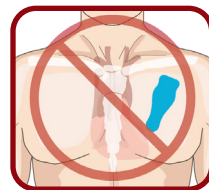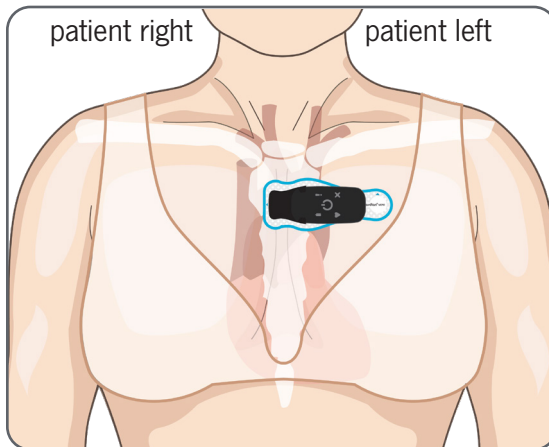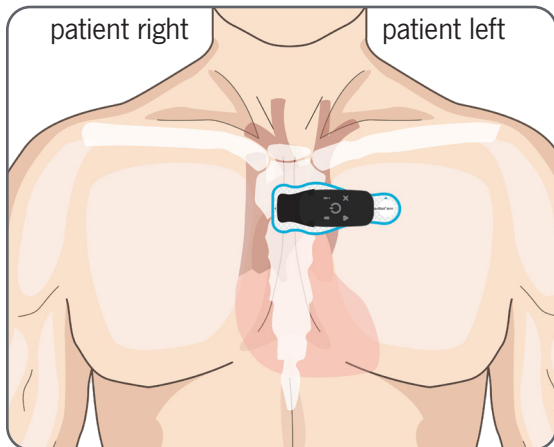

4

Turn on monitor

Press and release the center button  on the monitor once.

**If contact between the Strip and your skin is satisfactory:**

The monitor makes a quick series of three rising beeps and the green light blinks every 5 seconds. **You are successfully monitoring.**

**If contact between the Strip and your skin is NOT satisfactory:**

The monitor alerts you with two beeps every second. If this occurs:

1. Revisit the previous steps in this manual to ensure that:
 - The monitor is securely attached to the Strip (page 7).
 - BodyGuardian Mini is properly placed and securely adhered to your chest (pages 8-10).

When the monitor detects satisfactory skin contact, it will make a quick series of three rising beeps and the green light blinks every 5 seconds. **You are successfully monitoring.**

Did you receive Strips or an ECG lead set?

This section (pages 12-15) includes instructions for using an  ECG lead set to wear your BodyGuardian Mini monitor. If you received  Strips, see page 6.



Strip



ECG lead set

1 Prepare your skin

Familiarize yourself with the chest locations where you will attach ECG electrodes.

See **pages 14-15** for the chest locations where you will place the monitor and ECG electrodes.



If hair is present, remove the hair where you will place the monitor and ECG electrodes.






Use an enclosed saline wipe or soap and water to clean the area. Dry the area completely.



Do not apply lotions, oils, or perfumes.


Do not use alcohol wipes.

Do not abrade the skin.

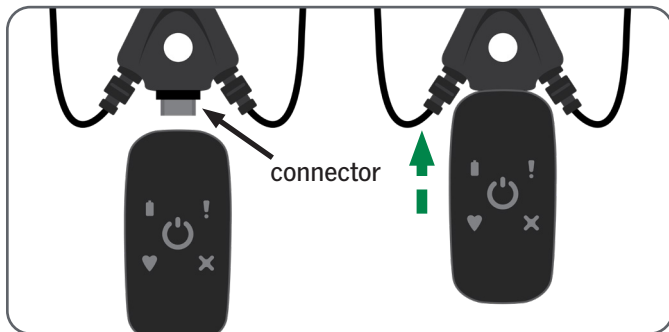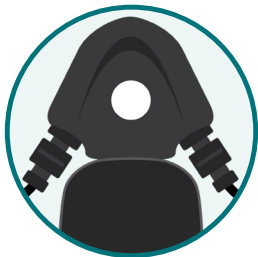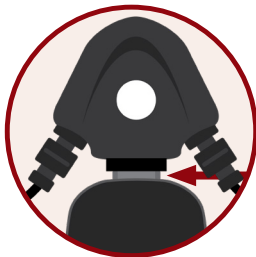**Good connection****Poor connection****2**

Attach the monitor and electrodes

1. Slide the BodyGuardian Mini monitor onto the center connector of the lead wires. Ensure the connector is inserted completely into the monitor.
2. Snap one ECG electrode into each of the leads, and one into the center monitor connector.

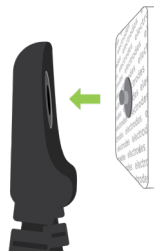




3

Attach to your chest

Remove the adhesive backing from each of the electrodes, then attach them to the following locations for the type of lead set you received with your equipment:

- 2-wire lead set
- 4-wire lead set (**next page**)


2-wire lead set

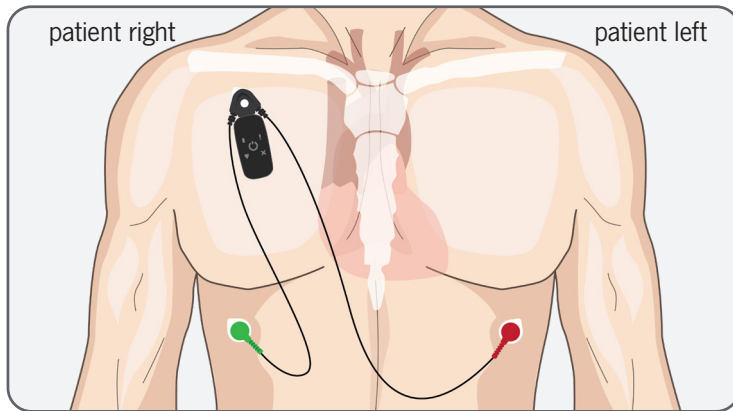

Apply only to intact, clean skin. Do not apply over open wounds, lesions, infected, or inflamed areas. If you have skin allergies, hypersensitivities, or experience skin irritation while wearing the monitor, consult your doctor to select the most appropriate option.



4-wire lead set

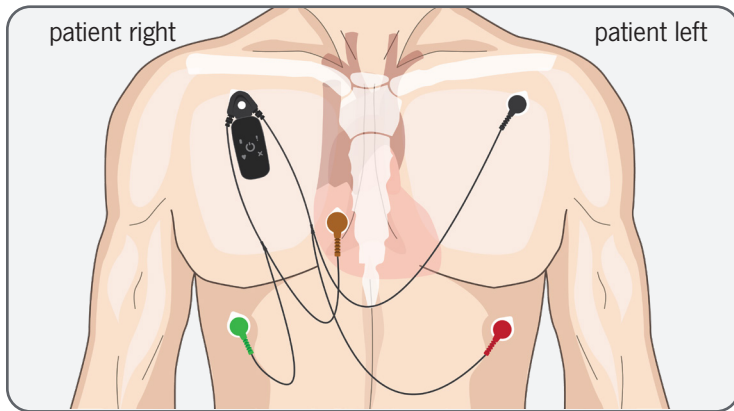

Apply only to intact, clean skin. Do not apply over open wounds, lesions, infected, or inflamed areas. If you have skin allergies, hypersensitivities, or experience skin irritation while wearing the monitor, consult your doctor to select the most appropriate option.



4

Turn on monitor

Press and release the center button  on the monitor once.

**If contact between the ECG lead set and your skin is satisfactory:**

The monitor makes a quick series of three rising beeps and the green light blinks every 5 seconds. **You are successfully monitoring.**

**If contact between the ECG lead set and your skin is NOT satisfactory:**

The monitor alerts you with two beeps every second. If this occurs:

1. Revisit the previous steps in this manual to ensure that:
 - The monitor and leads are securely attached to the ECG electrodes **(page 13)**.
 - The ECG lead set is properly placed and securely adhered to your chest **(pages 14-15)**.

When the monitor detects satisfactory skin contact, it will make a quick series of three rising beeps and the green light blinks every 5 seconds. **You are successfully monitoring.**

Recording your symptoms

If you begin to feel symptoms related to the reason your monitor was prescribed, such as **dizziness**, **chest pain**, or **shortness of breath**, you can manually record these events.

1. On the monitor, press and release the center button  once. The monitor makes a single beep on a successful button push. Do not hold the button down.
2. In your Holter diary, write down the symptom **date** and **time**.
3. Check the **symptoms** you are feeling.
4. Check the **activity** you were doing at the time you felt the symptom.

When you feel a symptom:

1. On the monitor, press and release the center button  once. Do not hold the button down.
2. Record the **DATE** and **TIME**.
3. Select at least one **SYMPTOM**.
4. Select your **ACTIVITY** at the time you felt the symptom.

	DATE	TIME	SYMPTOM (one or more)	None or accidental push	Light-headedness	Rapid or fast heartbeat	Flutter or skipped beats	Shortness of breath	Chest pain or pressure	Dizziness	Tired or fatigued	Passey
1	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
2	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
3	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
4	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
5	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
6	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
7	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐
8	___/___/___	___:___ AM/PM	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐



Be sure to include your diary in the box when you return your equipment at the conclusion of your prescription.

Charging the monitor



When the monitor battery is low, the monitor alerts you with two beeps and blue blinks every 5 seconds. If the low battery alert sounds, charge the monitor immediately.



Charge time: 1.5-2.5 hours



While charging, the blue light on the monitor will blink.

Once fully charged, the blue light will remain solidly lit.

1. Detach the monitor from the Strip or lead set on your chest. You do not need to remove the Strip or lead set from your chest when you charge the monitor.
2. Firmly plug the charger into the monitor.

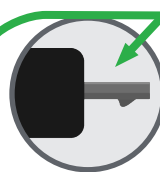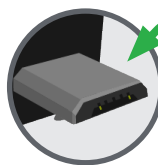

Important

Ensure the micro-USB connector is oriented this way when inserting into the monitor.



Do not use your own personal device chargers. Only use the charger(s) provided.

3. **AS SOON AS THE MONITOR IS FULLY CHARGED (1.5-2.5 HOURS)**, attach the charged monitor to the Strip or lead set on your chest.
4. Press and release the center button to power on the monitor.

If the monitor starts beeping during your monitoring session



2 beeps and 2 red flashes every 5 seconds

BodyGuardian Mini has not detected good contact with your skin for over 8 hours.

RESOLUTION

Revisit the previous steps in this manual to ensure that:

- The monitor is securely attached to the Strip (**page 7**) or ECG lead set (**page 13**).
- BodyGuardian Mini is properly placed and securely adhered to your chest (**pages 9-10, 14-15**).
- If the monitor continues to beep:
Remove the Strip or ECG electrodes, use a new Strip or ECG electrodes, and repeat the setup steps beginning on **page 6** (Strip) or **page 12** (ECG electrodes).



2 beeps and 2 blue flashes every 5 seconds

Your monitor battery is low.

RESOLUTION

Charge the monitor immediately (**page 18**).

Monitor indicators

* The monitor lights are faint. They can best be detected in a dark or dimly lit room.

Sounds		Lights*		What it means	Caused by
none			every 5 sec	Monitoring successfully.	
	and	none		Symptom marked.	Press-and-release button push.
	and		every 1 sec	Poor skin contact.	Monitor cannot detect good contact with skin to begin monitoring (p. 11, 16).
	and		every 5 sec	Poor skin contact.	Monitor has not detected good contact with skin for over 8 hours. Replace electrodes.
none	and		every 5 sec	Poor skin contact.	Monitor cannot detect good contact with skin during already established monitoring session. If this persists for 30 min or more, replace electrodes.
	and		every 5 sec	Monitor battery low.	Low monitor battery. Charge the monitor (p. 18).
	and			Monitor powering on.	Press-and-release button push.
	and	   		Monitor powering off.	Press-and-hold (9-12 sec) button push.
<i>After 9-12 seconds, all four monitor lights flash once, followed by a long beep.</i>					

Monitor indicators (continued)

Sounds	Lights*	What it means	Caused by
none and 	   	Device clock error.	Unable to start recording.
<i>Indicators are briefly lit, one at a time, in a rotating pattern.</i>			
none and    		Monitor has halted at internal error handling checkpoint.	Unable to start recording.
<i>Indicators blink.</i>			

Replacing the Strip

Change the Strip when it no longer adheres to the skin.

Strips are disposable in the trash. They are not biodegradable or compostable.



To protect signal quality, it's important to not expose the monitor's micro-USB connector directly to water.

1. Remove the monitor from the Strip, then charge the monitor.
2. Slowly peel back and roll the adhesive away from your chest, continually blotting with lukewarm soapy water.
TIP: For easiest removal, start by peeling back the edge of the Strip furthest from the black hub. It is suggested to remove the Strip while showering to limit the potential for skin irritation issues.
3. Use a *new* Strip to re-attach BodyGuardian Mini to your chest (see **pages 6-11**).

Replacing the ECG electrodes

Change the ECG electrodes every three days OR when they no longer adhere to the skin. Electrodes are disposable in the trash. They are not biodegradable or compostable.

1. Slowly peel back and roll the adhesive on the back of the electrodes away from your chest.
2. Remove the leads from the electrodes.
3. Remove the monitor from the center connector, then charge the monitor (**see page 18**).
4. Use *new* ECG electrodes to reattach BodyGuardian Mini to your chest (**see pages 12-16**).

**MRI**

BodyGuardian Mini cannot be worn while undergoing an MRI. Remove the monitor and the Strip or ECG lead set and electrodes. After MRI, re-apply the monitor to your chest.

**Air travel**

BodyGuardian Mini is safe to wear while traveling on aircraft.

Water resistance

The accessories (Strips and ECG lead set) are waterproof and can safely be submerged in water as deep as 3 feet. The monitor is waterproof when plugged into a body-worn accessory that does not leave the micro-USB port exposed. The micro-USB port should **not** be exposed to water.

- You do not need to remove BodyGuardian Mini when showering or bathing.
- If you will be swimming in water deeper than 3 feet, remove the monitor.
- If micro-USB connector on an unattached monitor is exposed to water, open air dry for at least 1 hour before reattaching to a Strip, lead set or charger.



After showering, bathing, swimming while wearing BodyGuardian Mini with a  Strip:

1. Pat the Strip with a towel to dry it off.
2. Press and hold the Strip or electrodes against your chest to ensure good contact with your skin.

RETURNING EQUIPMENT

Return your equipment immediately after your service ends


- ❗ With the device off your chest, turn off the monitor by pressing and holding the center button  for 9-12 seconds, until all 4 monitor lights illuminate and there is one long beep.

If returning your equipment to your health care provider:

1. Follow the instructions provided by your health care provider.

If returning your equipment to Boston Scientific Cardiac Diagnostic Services:

1. Place the powered-off monitor into the original postage-paid box. Then place the charging cord, Holter diary, unopened Strips, ECG lead set (if applicable), and patient survey into the box.
2. Remove the adhesive strip from flap on the shipping box and seal securely.
The prepaid return label is already affixed to the bottom of the box.
3. Drop package off at the shipping carrier listed on the prepaid return label.
For assistance, call Boston Scientific Cardiac Diagnostics at **888.500.3522**.

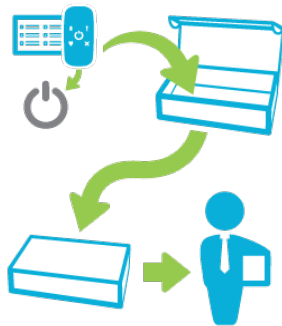

- ❗ If you have lost or damaged the equipment, Boston Scientific Cardiac Diagnostics will mail you an invoice for the cost.



On behalf of your physician, you are being provided with a heart monitor from Boston Scientific Cardiac Diagnostic Services. This test will be billed in two parts:

1. Your medical professional will bill your insurance for the in-office hook-up and the final reading of the test results.
 2. Boston Scientific Cardiac Diagnostic Services will bill your insurance for the use of the monitor along with providing the physician with 24/7 monitoring center service and all requested data and reports.
-



Your insurance company will usually send you an Explanation Of Benefits (EOB) describing the amount paid and the amount you owe Boston Scientific Cardiac Diagnostic Services. An EOB is not a bill.

Our goal is to help your physician diagnose your heart condition without delay or concerns about billing. We offer patient-friendly payment options, including a discounted rate at the time of service or a reduced rate when your insurance indicates a large fee payable by you for our services or the physician's services.



All equipment is the sole property of Boston Scientific Cardiac Diagnostic Services. To avoid financial liability and to ensure that you do not receive a bill for the value of the equipment in your possession (up to \$2,400), please return the equipment immediately after your service ends. Any request for an extension of your study beyond your prescribed study needs to come from your physician through a new order sent to Boston Scientific Cardiac Diagnostic Services. Be aware that the second study may not be covered by all insurance providers and you may be responsible for those charges.

Our Notice of Privacy Practices is available for your review at: **CDx.BostonScientific.com**



PATIENT SURVEY

We appreciate your feedback. Please return this survey inside the box once your study is completed. This survey can also be completed online at: **CDx.BostonScientific.com/us/en/patients/patient-support/patient-survey.html**

Patient name _____ City _____ State _____

Doctor's office / hospital _____ Contact me via: ☐ Phone _____

My hook-up took place at: ☐ Home ☐ Physician's office ☐ Email _____

BODYGUARDIAN™ MINI

	Exceeded expectations	Met expectations	Below expectations	Not applicable
Rate our service level	1	2	3	4
General satisfaction with Boston Scientific Cardiac Diagnostics	☐	☐	☐	☐
Ability to contact us at 888.500.3522	☐	☐	☐	☐
Explanations given to hook up your monitor	☐	☐	☐	☐
Amount of time the representative spent with you	☐	☐	☐	☐
Knowledge, professionalism, courtesy of phone representatives	☐	☐	☐	☐
Rate the contents of the materials	1	2	3	4
Printed instruction manual included in your box	☐	☐	☐	☐

continue ➔

PATIENT SURVEY

BODYGUARDIAN™ MINI

	Excellent	Average	Poor	Not applicable
Rate your experience using the monitor and accessories	1	2	3	4
Setting up the monitor	☐	☐	☐	☐
Wearing and charging the monitor	☐	☐	☐	☐
Adherence and comfort of the Strips	☐	☐	☐	☐
Shipment of supplies for the length of your study	☐	☐	☐	☐

Comments

☐ I am interested in potentially sharing my patient experience and consent to Boston Scientific Cardiac Diagnostics, Inc. contacting me for such purposes at the contact information I provided above, which may include telephone or email.

